

ATMOS VideoScope

USB Videoendoscope for Swallow Studies

*Exactly
what
I need*



ATMOS VideoScope


Cutting edge technology and perfect handling—plug and play.

The **ATMOS VideoScope** combines all important features for stationary and mobile examinations. Portability and handling are unique: the light source and camera processor are integrated in the endoscope itself. To start your FEES study, just check that the USB endoscope is connected to your PC and then launch the video software. Crisp images are produced thanks to distal technology, and two diameter options down to 2.9 mm allow gentle examinations even in pediatrics.

In addition to the integrated camera processor, the **ATMOS VideoScope** endoscope also features an integrated high-power LED that produces the brightness required for clear vision and accurate diagnosis.



Simple and safe reprocessing is an essential aspect of FEES examinations. Manual cleaning procedures are both gentle and safe. The **ATMOS VideoScope** can be completely immersed into liquids.

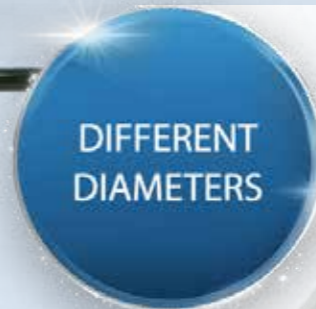

The **ATMOS VideoScope** is the result of highly innovative technology. Its various diameters makes it the perfect choice for scoping both children and adults.



Distal chip technology stands for high-resolution image and small endoscope diameter at the same time. The **ATMOS VideoScope** produces crystal clear images with a significantly higher resolution than that of standard fiberoptic endoscopes.





The **ATMOS Capture Suite** software allows you to record and store videos as well as to view them later and to fulfill your reporting duties. Still image and video recording can be initiated via the capture button on the endoscope.

Technical Data

	VideoScope	VideoScope PED
Diameter	3.9 mm	2.9 mm (PED)
Working length		310 mm
Field of vision		90°
Angulation		130°/130°
Sensor resolution		400 x 400 pixel
Weight		330 g
Interface		USB
Light source		Integrated LED

Wiring Diagram

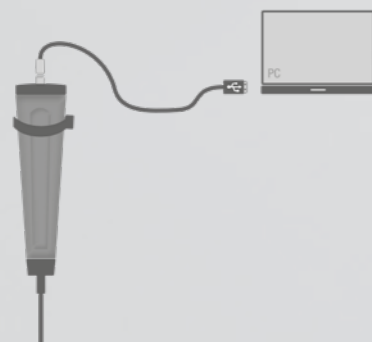

Quality  Made in Germany





MedizinTechnik

ATMOS MedizinTechnik GmbH & Co. KG

Ludwig-Kegel-Str. 16
79853 Lenzkirch
Germany

Tel +49 (0)7653-689-0
Fax +49 (0)7653-689-190
E-Mail info@atmosmed.de

www.atmosmed.com