



---

ATMOS C 11 Systema

# The treatment unit for work efficiency

---

*Exactly  
what  
I need*

---

# C 11 Systema

The little giant—small for starters, but ready to grow

Every ENT practice is different. Thanks to its intelligent modular system, the ATMOS C 11 Systema is the right solution for any room situation and budget. Its compact design makes the most out of even the smallest spaces. High-quality materials as well as sophisticated workspace organization were taken into account when developing the ATMOS C 11 Systema.



## Maximum Flexibility

### THANKS TO COMPREHENSIVE CONFIGURATION OPTIONS

- Modular configuration
- Full range of functions even for the smallest of treatment units

## Work efficiency

### THANKS TO AN INNOVATIVE OPERATING CONCEPT

- Control element positions that optimize your workflow
- Autostart function

## Clear overview

### THANKS TO THE NEAT AND STRUCTURED WORKSPACE

- Instrument management on the work surface
- Drawers for storage space
- Various holders

info



Which treatment unit is the right choice for my requirements?

ATMOS offers ENT doctors different treatment units that can transform their individual ideas into reality. You can find more information on our treatment units on our website.

Designed,  
developed and  
made in Germany



ATMOS MODULAR SYSTEM

# Configuration? Comprehensive.

The ATMOS C 11 Systema provides you with all the important elements for ENT treatments in one compact treatment unit. The sky's the limit when it comes to selecting the modules and devices you need for carrying out your examinations and treatments. Configure the treatment unit with our modular system to suit your individual requirements, space, and processes.

### LIGHT MODULES

Our LED light modules are cooled passively by a sophisticated ventilation system, are completely noiseless, and don't stir up any dust. Thanks to the optimized light connection system, the endoscope and tissue never heat up, even when the module is running at maximum power.



Scan the QR code for detailed information

### VISUALIZATION SYSTEMS

Endoscopy and visualization make ENT diagnostics more efficient and help patients to understand their diagnoses. You, too, can benefit from our innovative imaging diagnostics solutions. Thanks to the various connection possibilities of our devices, they can be integrated easily. Our ATMOS Capture Suite visualization and archiving software provides additional process optimization.



### DIAGNOSTIC EQUIPMENT

Optimize your workflow for routine examinations and cut your patients' waiting times with our ENT diagnostics package.



## Maximum flexibility

**BOTTLE MANAGEMENT**  
Holders for ear irrigation bottles and sprayer bottles

**ENDOSCOPE QUIVERS**  
Ergonomic and hygienic holders

**LED POWER SOURCES**  
LED connections for headlamp and LED light handle for direct connection to optical devices



**INSTRUMENT MANAGEMENT**  
All instruments are within immediate and easy reach

**DRAWERS**  
For storing consumables and rarely needed instruments

**HANDLE AND HOSES**  
Holders for suction hose and compressed air handle



**METAL BODY**  
Corrosion-resistant metal body with easy-to-clean surfaces

**SUCTION PUMP**  
Integrated suction module with autostart function



**SECRETION MANAGEMENT AND DDS**  
Easy manual emptying of secretion canisters by means of the ATMOS Direct Docking System (DDS)



ATMOS WORKFLOW

# Sustainable concept. Many possibilities.

A lack of space in outpatient clinics or small medical practices often means that not all the necessary functions can be placed within reach and in an ergonomic position. With the ATMOS C 11 Systema, however, a treatment unit awaits you that, thanks to its compact design, satisfies both ergonomic and process-optimized needs and fits perfectly into your ENT environment.

Work efficiency



### AUTOMATIC ACTIVATION

Enjoy a more relaxed and efficient way of working: all functions of the treatment unit are started automatically when the handles are taken out of their holders and can be operated with one hand. The optimum positioning of the control elements also allows you to optimize your daily workflow.



### DIFFERENT WORKING HEIGHTS

The ATMOS C 11 Systema is available in two versions: sitting and standing. The ATMOS C 11 Systema (sitting) is intended for treatment of patients while the doctor is sitting and therefore has a lower working height. The ATMOS C 11 Systema (standing), on the other hand, is designed for doctors who treat their patients while standing.

info

### ATMOS CHAIR COMFORT ATMOS CHAIR PROFESSIONAL

The central element in an ENT examination room is the patient chair: make sure your patients feel safe and comfortable with an ATMOS Chair.

A modern design and a great feel thanks to the soft and extremely durable upholstery material are the new standard for all ATMOS Chairs.



Scan the QR code to go to the ATMOS Treatment Chairs web page.



## ATMOS INSTRUMENT MANAGEMENT

# Comprehensive solutions. Maximum process safety.

When it comes to instrument management, benefit from the clever solutions and special features for expanding the ATMOS C 11 Systema. In addition to the corrosion-resistant metal cabinet and easy-to-clean surfaces, it is in particular the large work surface, drawers, and various storage possibilities that enable hygienic working. The flexibility in selecting only the modules you need for your everyday work leaves nothing to be desired, also with regard to instrument management.

## INSTRUMENT MANAGEMENT



With the ATMOS C 11 Systema, instrument management is located directly on the work surface. All important instruments are therefore within immediate and easy reach. Four large trays (284 x 184 mm) are optionally available on the upper level for instrument deposit.

**Enhanced features:** By expanding the instrument deposit to two levels, you can create space for a total of six large trays. The lower level can be pulled out to open it.

Clear  
overview

## DRAWERS



One small drawer and several large drawers for storing consumables and instruments that are rarely used provide for a neat and structured workspace.

**Enhanced features:** Optionally, you can equip your drawer with an instrument disposal. When the drawer is closed, the instruments can be dropped into the integrated plastic container (instrument disposal tray). The drawer can be opened to safely replace the instrument disposal tray.

## MEDICATION SPRAYERS AND EAR IRRIGATION



Holders for medication sprayers, ear irrigation bottles, and endoscope quivers are integrated directly into the work surface.

The additional compressed air module enables you to apply various types of medications (aqueous, oily, powdery) by connecting the sprayer bottles. In combination with the ear irrigation module, warmed water can be used for ear irrigation using compressed air. All water-bearing parts are autoclavable and therefore particularly easy to reprocess. Smart endoscope management allows you to easily store clean and used endoscopes separately.

## SUCTION AND SECRETION MANAGEMENT



The quiet suction system has a user-friendly autostart function that starts the system when the suction hose is removed from the holder. Furnished with a suction module, the ATMOS C 11 Systema is optimally prepared for all ENT applications requiring suction.

Thanks to the easy emptying of the secretion canister by means of the ATMOS Direct Docking System, suction becomes even more user-friendly. The secretion canister is equipped with a mechanical overflow protection and a hydrophobic bacterial filter. As an alternative, various disposable canister systems can also be used.

## ATMOS C 11 SYSTEMA

# Treatment unit? ATMOS C 11 Systema.


The ATMOS C 11 Systema contains exactly what ENT doctors need to effectively and efficiently meet the various demands arising in highly frequented offices: nothing more, but also nothing less. Its compact design makes the most out of even the smallest spaces.



## ATMOS C 11 SYSTEMA — TECHNICAL DATA



	ATMOS C 11 Systema (sitting)	ATMOS C 11 Systema (standing)
Dimensions (H x W x D)	840 x 960 x 535 mm	960 x 960 x 535 mm
Large drawer (600 x 120 mm)	2x	3x
Small drawer (314 x 148 mm)		1x
Instrument disposal		optional
Side pull-out for secretion canister		✓
Suction		✓
LED power sources		2x
Holders for handles and hoses		✓
Holders for endoscope quivers		✓
Holders for sprayer bottles		✓
REF	541.0000.0	542.0000.0

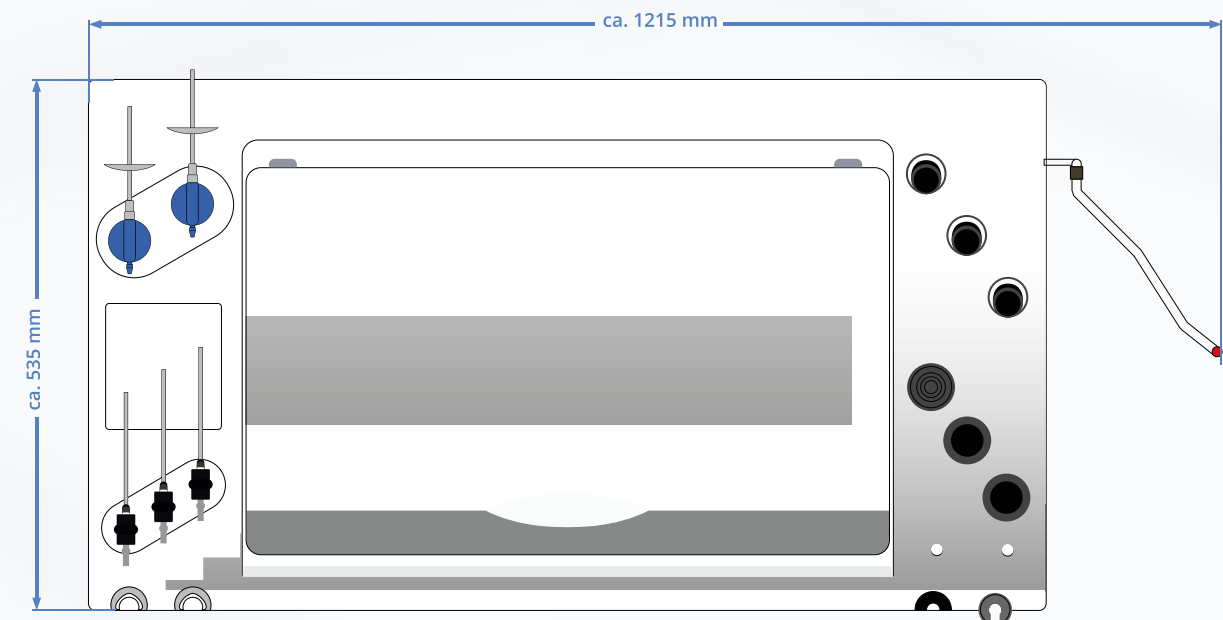



**ATMOS MedizinTechnik GmbH & Co. KG**

Ludwig-Kegel-Str. 16  
79853 Lenzkirch / Germany

Phone +49 (0)7653-689-0  
Fax +49 (0)7653-689-190  
E-mail [info@atmosmed.com](mailto:info@atmosmed.com)

**[www.atmosmed.com](http://www.atmosmed.com)**



Scan the QR code  
to go to the  
ATMOS C 11 Systema  
web page.